

# ECO ASSIST

EMPOWER YOUR MOBILITY



REHASENSE



Exclusively distributed by

**aidacare™**  
HEALTHCARE EQUIPMENT


# ECO ASSIST

## LIGHTWEIGHT POWER ADD-ON

ECO ASSIST is a powered add-on that transforms a manual wheelchair into a motorised mobility aid. Made from aluminium and steel, its lightweight and durable design is easy for users to dock to their wheelchair in seconds. Featuring a robust motor powered by an efficient lithium-ion battery, it's ideal for navigating city streets with ease. With adjustable controls and reliable braking, ECO ASSIST allows users to enjoy everyday journeys with peace of mind.



**Dual Brake System**  
Mechanical disc and electric braking system

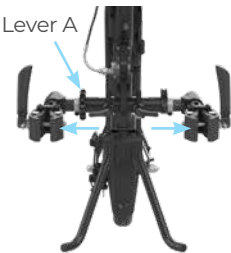

**Width Adjustable**  
Quick adjustment via manual clamps



**Leg Supports**  
Height adjustable between 260-420mm



**Quick Disconnect**  
Supports easier transportation and storage with a single point disconnection

- Easy manual docking system allows for quick and easy mounting
- Compatible with both rigid and folding wheelchairs, with fixed front frames or removable legrests
- Universal clamp system allows fitting to most wheelchair types\*
- 160mm disc brake for strong, reliable braking
- Fits wheelchairs with a width of 210 – 480mm (between footplates)
- Clamp height adjustable from 260 – 420mm
- Walk-assist mode
- Traction control for improved grip
- Reverse gear for easier manoeuvring

### Accessories



**Additional Footrest**  
WCA872620



**Bag**  
WCA872268



**2 YEAR WARRANTY**

**Airline Approved Battery Pack - 5.6AH**  
WCA872235

Eco Assist	Wheel Size	Engine Power	Max. Engine Torque	Speed Modes	Max Speed	Driving Range	Overall Weight	Weight of Heaviest Part	SWL
WCA872000	12"	350W	39Nm	3	10km/h	Up to 30km	15.5kg	10kg	110kg

\*Not to be used with Carbon Fibre wheelchair frames; warranty will be voided.